




**KOALA**  
Medical Pty Ltd

EST. 1993

# VentiSeal Endotracheal Tubes

Efficient, low pressure seal, available in a wide range of shapes and sizes for all airway needs.

A thin ultra smooth High Volume Low Pressure Cuff, that effectively conforms to the tracheal wall.

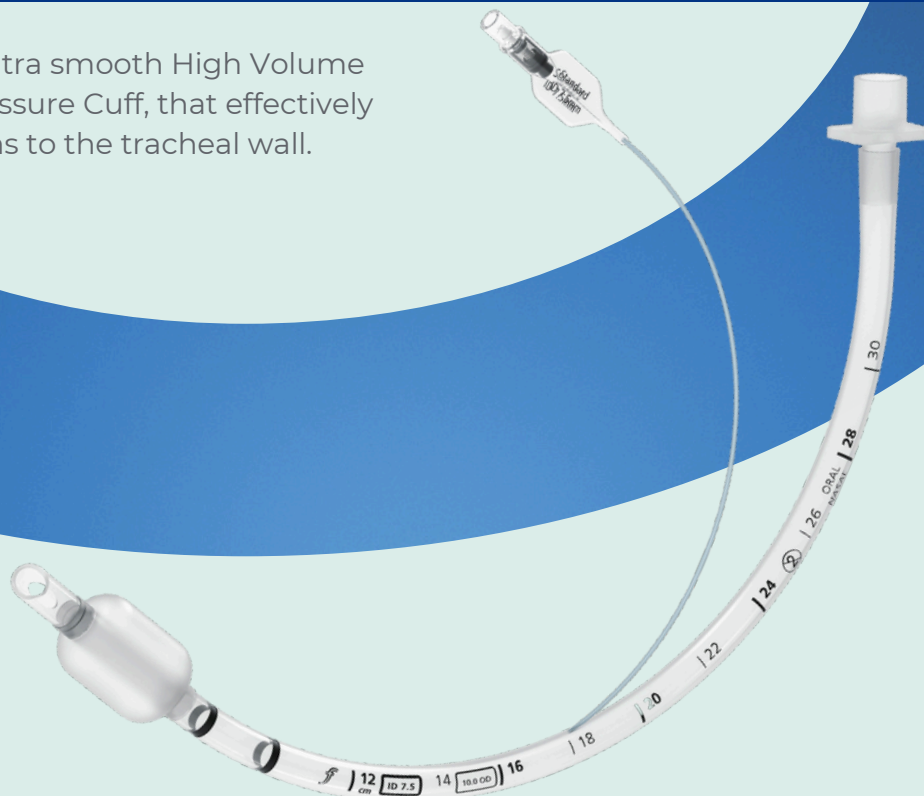

- + Sterile
- + Single use
- + Latex free
- + DEHP free



## Soft Rounded Distal Tip

Gentle, bevelled tip reduces potential for soft tissue trauma during insertion, with minimal increase in flow resistance.



## Large Murphy's Eye

Allows larger flow of gas, to occur, should the tip become occluded.



## Circumferential Vocal Cord Markers

Double guide lines are clearly seen from all angles, facilitating safe and accurate cuff placement below the vocal cords.



02 9979 9374



[sales@koalamedical.com.au](mailto:sales@koalamedical.com.au)



[koalamedical.com.au](http://koalamedical.com.au)





**KOALA**  
Medical Pty Ltd

EST. 1993

# VentiSeal Endotracheal Tubes

Efficient, low pressure seal, available in a wide range of shapes and sizes for all airway needs.



## Distance Markings

Clear and bold print give earlier indication of depth and is a simple, safe and reliable method for correct tracheal tube placement.



## Kink Resistant

Retains tube patency for patient safety, whilst softening at body temperature to conform to patient's respiratory tract.



## Coloured Inflation Line

For quick identification to differentiate from other tubing.



## Simple Verification

Small rounded shape of the balloon gives a more responsive and tactile feel to the user for inflation check.

## Reinforced Central Rib

Internal line running through the centre of the pilot balloon prevents complete collapse during rapid deflation. Line to cuff is always open.

ORAL/NASAL	CURVED (ORAL) SOUTH	CURVED (NASAL) NORTH	DESCRIPTION	QUANTITY
FL-038-971-025			Endo tube 2.5mm	10 pack
FL-038-971-030	FL-038-972-030	FL-038-974-030	Endo tube 3.0mm	10 pack
FL-038-971-035	FL-038-972-035	FL-038-974-035	Endo tube 3.5mm	10 pack
FL-038-971-040	FL-038-972-040	FL-038-974-040	Endo tube 4.0mm	10 pack
FL-038-971-045	FL-038-972-045	FL-038-974-045	Endo tube 4.5mm	10 pack
FL-038-971-050	FL-038-972-050	FL-038-974-050	Endo tube 5.0mm	10 pack
FL-038-971-055	FL-038-972-055	FL-038-974-055	Endo tube 5.5mm	10 pack
FL-038-971-060	FL-038-972-060	FL-038-974-060	Endo tube 6.0mm	10 pack
FL-038-971-065	FL-038-972-065	FL-038-974-065	Endo tube 6.5mm	10 pack
FL-038-971-070	FL-038-972-070	FL-038-974-070	Endo tube 7.0mm	10 pack
FL-038-971-075			Endo tube 7.5mm	10 pack
FL-038-971-080			Endo tube 8.0mm	10 pack
FL-038-971-085			Endo tube 8.5mm	10 pack
FL-038-971-090			Endo tube 9.0mm	10 pack



02 9979 9374



sales@koalamedical.com.au



koalamedical.com.au

